


# Metal-free Finishing and Polishing with Flexible Disks

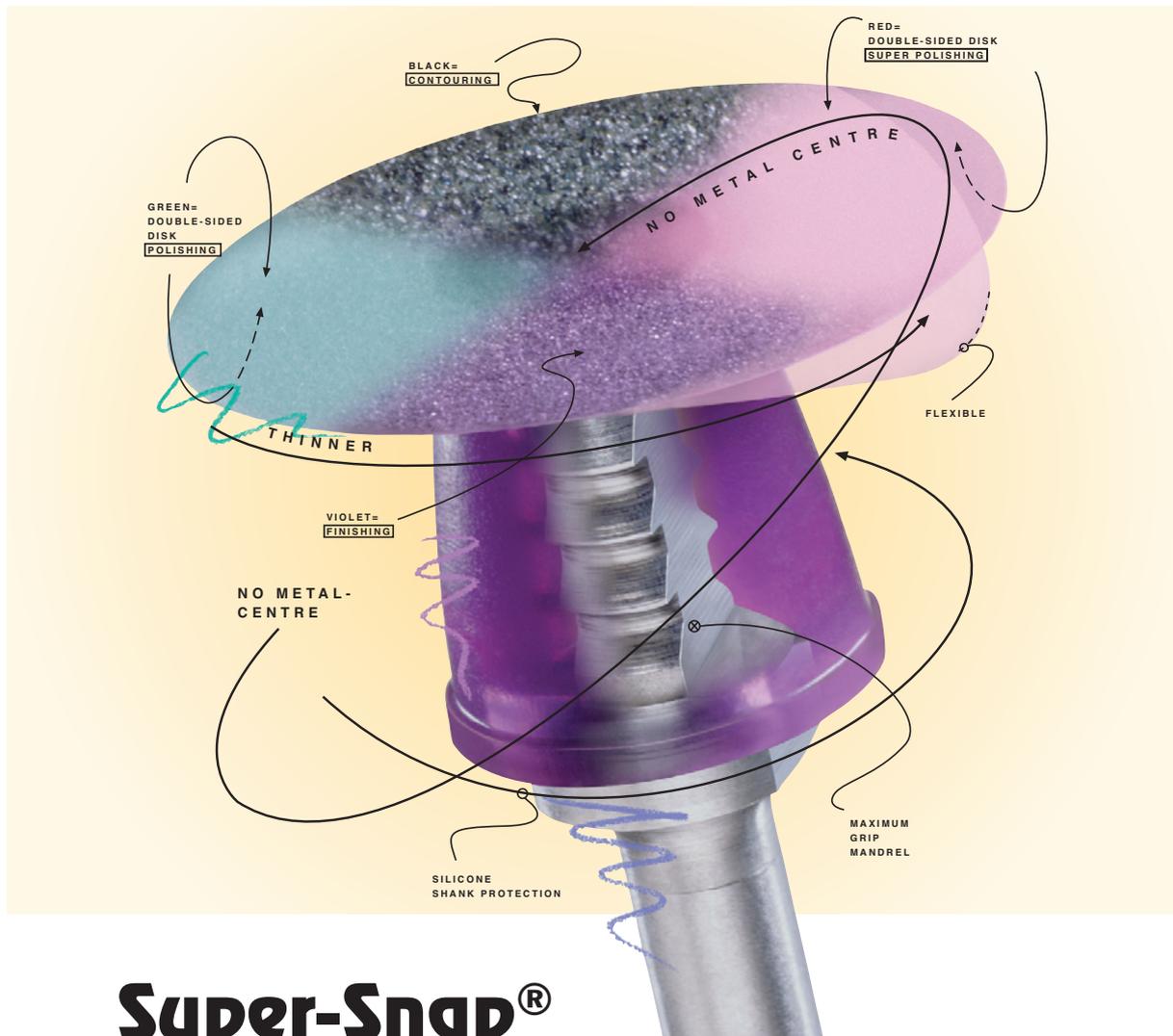

## Super-Snap®

**Super-Snap** disks are designed for easy, fast and safe contouring, finishing and polishing of all microfilled and hybrid composites. The elastic shank mount allows easy placing of the disk onto the mandrel.

Due to the fact that the disks have no metal centre, gouging and discolouration of the composite is avoided. The ultra-thin disks are very flexible, bending under pressure without breaking and allowing easier access to interproximal areas.

Colour-coordinated in four different grits, the complete range of disks can be found in the **Super-Snap** Rainbow Technique Kit.




Proven Products for Better Dentistry

# Step-by-Step to the Finest Restorations

- **Colour-coordinated**  
Clear, logical and ergonomic arrangement of the instruments
- **Black – Coarse Disk**  
Uniform coarse grit particles for fast, smooth composite contouring
- **Thin, Double-sided Disks**  
For safe, fast polishing without changing or inverting disks
- **Secure Mandrel Design**  
Holds disks firmly for maximum control during polishing
- **No Metal Centre on the Polishing Surface**  
Prevents gouging or discolouration of the composite surface
- **Thin Polystrips**  
For easier access to the approximal areas; these are colour-coordinated to the Super-Snap disks with the same grits – black/violet (coarse and medium), green/red (fine and superfine)


## Surface Roughness (Ra)

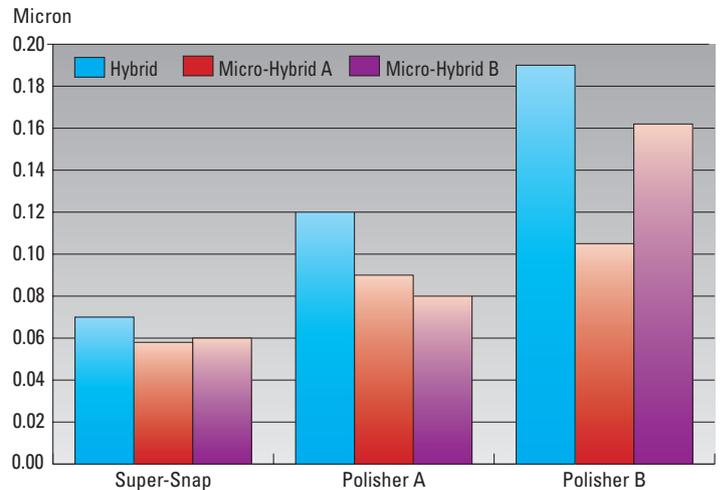

Compared to other polishing disks, the **Super-Snap** system achieves excellent results especially on hard-to-polish composites.

## Super-Snap Details and Order Information:

	Contouring			Finishing			Polishing		Superpolishing	
<b>Super-Snap Disks</b>										
(50 pcs./box)	<b>Black</b> (coarse)			<b>Violet</b> (medium)			<b>Green</b> (fine)		<b>Red</b> (superfine)	
<b>Super-Snap Mandrel</b>										
	FG – 0438 (6 pcs./box)			CA – 0439 (6 pcs./box)			CA – 0440 Plastic (100 pcs./box)			
<b>Super-Snap Polystrips</b>										
(100 pcs./box)	<b>Black/Violet</b> (coarse/medium)						<b>Green/Red</b> (fine/superfine)			

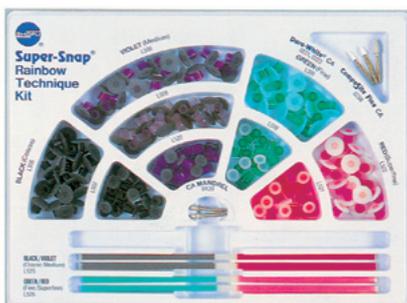

PN 0500

The **Super-Snap Rainbow Technique Kit** contains 180 disks (100 pcs. standard, 80 pcs. mini), 40 polystrips, 2 **Dura-White** stones (CN1 and FL2), 1 **CompoSite** Fine midi-point, 4 steel mandrels, CA shank



## SHOFU DENTAL GMBH

Am Brüll 17 · 40878 Ratingen / Germany  
Phone: +49 (0) 21 02 / 86 64-0 · Fax: 86 64-64  
E-mail: info@shofu.de · www.shofu.de